



FILTER TECHNOLOGY

CARTFLOW

Junior Cartridge

CATALOG





INDEX

**General Application Junior Cartridge filters
PES membrane**

2

**Bio-Burden Reduction Junior Cartridge filters
PES membrane**

4

**Sterilizing Junior Cartridge filters
PES membrane**

6

**Bio-burden Reduction Junior Cartridge filters
Hydrophobic PTFE membrane**

8

**General Application Junior Cartridge filters
Hydrophobic PTFE membrane**

10

**General Application Junior Cartridge filters
Nylon6,6 membrane**

12


**General Application Junior Cartridge filters
Hydrophilic PTFE membrane**

14

**Junior Cartridge filters
Polypropylene media**

16

Drawings of Junior Cartridge Filters

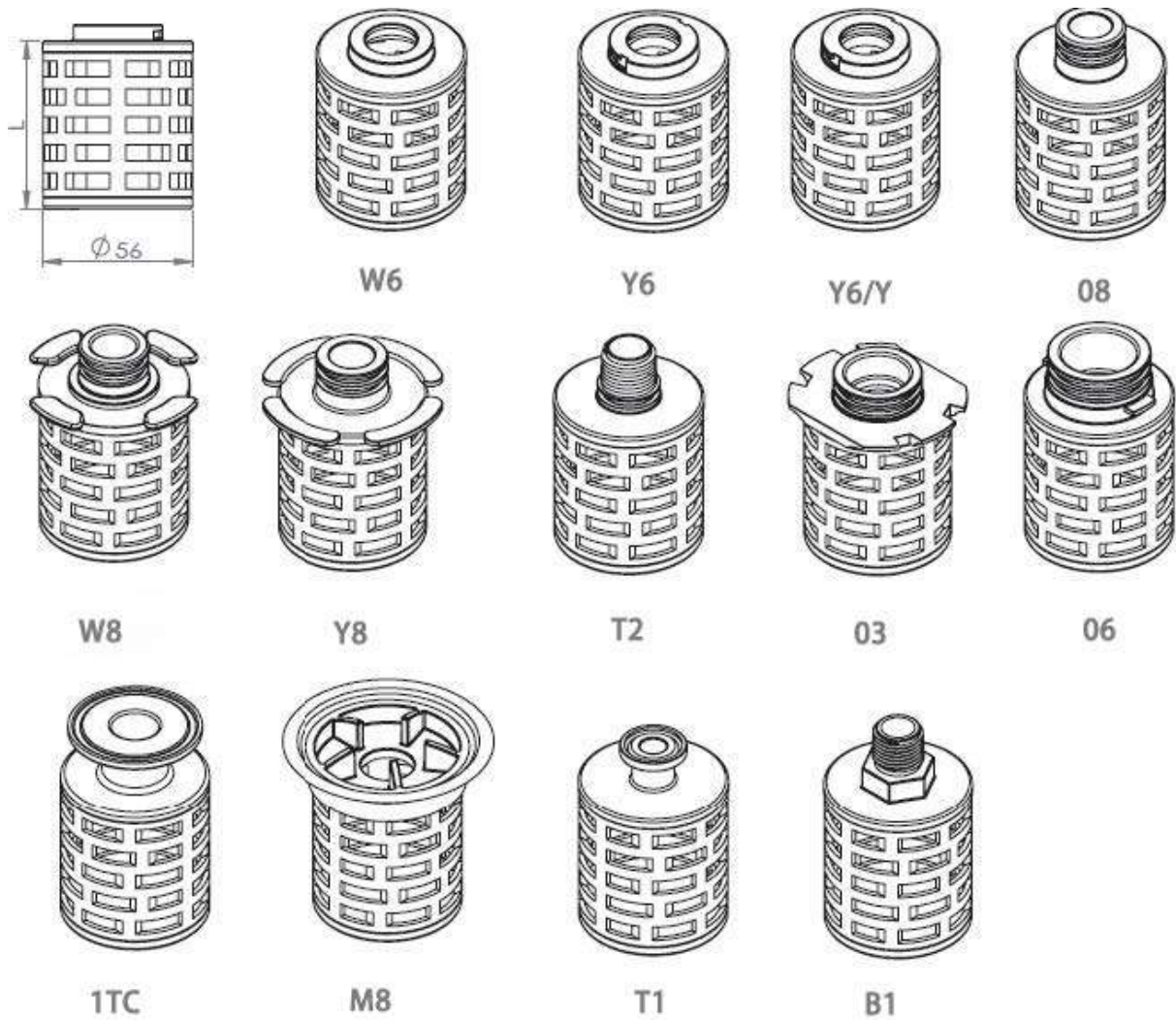

L: available in

32mm,63mm,
68mm,94mm,124mm

Fittings:

W6=116 w/o locking groove
 Y6=116 w/ locking groove
 Y6/Y=116 w/ locking groove with SS reinforced
 08=118 double o-ring
 W8=118 w/ short lugs
 Y8=118 w/longlugs
 T2=1/2"NPT threaded
 03=123 double o-ring
 06=126 double o-ring
 1TC=1"TC
 M8= 80mm skirting
 T1=1/2"TC
 B1=1/2"BSP-T

General Application Junior Cartridge filters/PES membrane

CMK sterilizing Grade Filter Polyethersulfone Cartridge are 100% integrity tested, it is constructed of single layer hydrophilic asymmetric polyethersulfone membrane. It offers broad chemical compatibility, higher filter area with high flow rates at low pressure drops, and low extractables.



Typical Applications

- Bottled Water
- Process water
- Vinegar
- Aseptic Packaged Liquids
- Beer, Wine and Spirits
- Juice, Soft drinks, Edible Oils
- Optical Disk manufacturing
- Hard dDsk manufacturing
- Optical Coatings
- Bulk Chemicals
- Pharmaceutical intermediates

Filter Area

HH: 500cm² SS: 1000cm²
DD: 1500cm² LL: 2100cm²

Integrity Test Specifications

Water wetted membrane

Pore size	Minimum Bubble point
0.1um	70 psig (4.8 bar)
0.2um	45 psig (3.1 bar)
0.45um	24 psig (1.7 bar)
0.65um	19 psig (1.3 bar)
0.8um	17 psig (1.2 bar)
1.2um	11.6 psig (0.8 bar)

Sanitization/Sterilization

- Filtered Hot water 90°C(194°F), 30 minutes, multiple cycles, max 3 psid forward flow
- Autoclave 121°C(250°F), 30 min, multiple cycles
- In-line Steam 135°C(275°F), 30 min, multiple cycles

Maximum Operating Parameters

- Differential Pressure
 - Forward 50 psid (3.4 barg) at 20°C (68°F)
 - Reverse 40 psid (2.7 barg) at 20°C (68°F)
- Operating Temperature 82°C(180°F) at 10 psid (0.69 barg) in water
- Recommended Changeout Pressure: 35psid(2.4barg)

Toxicity

All components meet the specifications for biological safety per USP Class VI-121°C for plastics.

Construction Materials

Filter Membrane: Polyethersulfone (PES) Membrane

Media Support: Polypropylene

End Caps: Polypropylene

Inner Core: Polypropylene

Outer Cage: Polypropylene

Sealing Method: Thermal Bonding

O-rings: Buna, Viton(or FKM), EP, Silicone,

Encapsulated Silicone, Encapsulated Viton (or FKM)

Food Safety Compliance

Material of construction comply with FDA regulations for food and beverage contact use as detailed in the US Code of Federal Regulation, 21CFR.

Material used to produce filter media and hardware are safe for use in contact with foodstuffs in accordance with EU Directives 10/2011.

Eg.=> CMKPS0010GHHOSW6PBY0

ORDERING INFORMATION										
Product Type	Membrane Type	Removal Rating	Application	Size	Diameter	Endcap	Inner Core	Sealing Material	Connection Support	Revision
CMK= Mini Pleated Cartridge	PS = PES	0010=0.1um	G = Gen Purpose	HH=500cm ²	OS=STD	W6=116 w/o locking groove Y6=116 w/ locking groove O8=118 double o-ring W8=118 w/ short lugs Y8=118 w/longlugs T2=1/2"NPT threaded O3=123 double o-ring O6=126 double o-ring 1TC=1"TC M8= 80mm skirting T1=1/2"TC B1=1/2"BSP-T	P=PP	B=Buna-N	Y= With SS reinforced	0 = Rev.0
		0020=0.2um		SS=1000cm ²				E=EPDM		
		0045=0.45um		DD=1500cm ²				S=Silicone		
		0065=0.65um		LL=2100cm ²				V=Viton		
		0080=0.8um								
0120=1.2um										

Bio-Burden Reduction Junior Cartridge filters/PES membrane

CMK Bio-Burden Reduction Grade Filter Polyethersulfone Junior Cartridge are 100% integrity tested; ensures bio-burden and small particle removal in a wide range of food&beverage, biological liquids, and intermediate bulk pharmaceutical fluids. The CMK Series is constructed of single layer hydrophilic asymmetric polyethersulfone membrane. It offers broad chemical compatibility, high flow rates at low pressure drops, and low extractables.



Typical Applications

- Cell Culture Media
- Large Volume Parenterals(LVP's)
- Pharmaceutical Bulk Chemical Solutions
- Diagnostics
- Blood and Serum Fractions
- Purified Water
- Beer, Wine and Spirits
- Juice & Soft Drinks
- Bottled Water

Filter Area

HH: 500cm² SS: 1000cm²
DD: 1500cm² LL: 2100cm²

Integrity Test Specifications

Water wetted membrane

Pore size	Minimum Bubble point
0.2um	3.5barg
0.45um	2.3barg
0.65um	1.5barg

Maximum Operating Parameters

- Differential Pressure
 - Forward 50 psid (3.4 barg) at 20°C (68°F)
 - Reverse 40 psid (2.7 barg) at 20°C (68°F)
- Operating Temperature 82°C(180°F) at 10 psid (0.69 barg) in water
- Recommended Changeout Pressure: 35psid(2.4barg)

Sanitization/Sterilization

- Filtered Hot water 90°C(194°F), 30 minutes, multiple cycles, max 3 psid forward flow
- Autoclave 121°C(250°F), 30 min, multiple cycles
- In-line Steam 135°C(275°F), 30 min, multiple cycles

Toxicity

All components meet the specifications for biological safety per USP Class VI-121°C for plastics.

Construction Materials

Filter Membrane: Polyethersulfone (PES) Membrane

Media Support: Polypropylene

End Caps: Polypropylene

Inner Core: Polypropylene

Outer Cage: Polypropylene

Sealing Method: Thermal Bonding

O-rings: Buna, Viton(or FKM), EP, Silicone,
Encapsulated Silicone,
Encapsulated Viton (or FKM)

Food Safety Compliance

Material of construction comply with FDA regulations for food and beverage contact use as detailed in the US Code of Federal Regulation, 21CFR.

Material used to produce filter media and hardware are safe for use in contact with foodstuffs in accordance with EU Directives 10/2011.

Eg.-> CMKPS0020BHHOSW6PBY0

ORDERING INFORMATION										
Product Type	Membrane Type	Removal Rating	Application	Size	Diameter	Endcap	Inner Core	Sealing Material	Connection Support	Revision
CMK=	PS = PES	0020=0.2um	B =	HH=500cm ²	OS=STD	W6=116 w/o locking groove	P=PP	B=Buna-N	Y= With SS	0 = Rev.0
Mini		0045=0.45um	Low	SS=1000cm ²		Y6=116 w/ locking groove		E=EPDM	reinforced	
Pleated			Bio			O8=118 double o-ring				
Cartridge		0065=0.65um		DD=1500cm ²		W8=118 w/ short lugs		S=Silicone		
				LL=2100cm ²		Y8=118 w/longlugs		V=Viton		
						T2=1/2"NPT threaded				
						O3=123 double o-ring				
						O6=126 double o-ring				
						1TC=1"TC				
						M8= 80mm skirting				
						T1=1/2"TC				
						B1=1/2"BSP-T				

Sterilizing Junior Cartridge filters/PES membrane

CMK sterilizing Grade Filter Polyethersulfone Cartridge are 100% integrity tested; it is constructed of double layer hydrophilic asymmetric polyethersulfone membrane. It offers broad chemical compatibility, high flow rates at low pressure drops, and low extractables.



Typical Applications

- Cell Culture Media
- Large Volume Parenterals(LVP's)
- Pharmaceutical Bulk Chemical Solutions
- Diagnostics
- Blood and Serum Fractions
- Purified Water
- Beer, Wine and Spirits
- Juice & Soft Drinks
- Bottled Water

Filter Area

HH: 500cm² SS: 1000cm²
DD: 1500cm² LL: 2100cm²

Integrity Test Specifications

Water wetted membrane

Pore size	Minimum Bubble point
0.1um	2.2barg(water/IPA)
0.2um	3.8barg(water)
0.45um	2.7barg(water)

Maximum Operating Parameters

- Differential Pressure
 - Forward 50 psid (3.4 barg) at 20°C (68°F)
 - Reverse 40 psid (2.7 barg) at 20°C (68°F)
- Operating Temperature 82°C(180°F) at 10 psid (0.69 barg) in water
- Recommended Changeout Pressure: 35psid(2.4barg)

Sanitization/Sterilization

- Filtered Hot water 90°C(194°F), 30 minutes, multiple cycles, max 3 psid forward flow
- Autoclave 121°C(250°F), 30 min, multiple cycles
- In-line Steam 135°C(275°F), 30 min, multiple cycles

Toxicity

All components meet the specifications for biological safety per USP Class VI-121°C for plastics.

Construction Materials

Filter Membrane: Polyethersulfone (PES) Membrane

Media Support: Polypropylene

End Caps: Polypropylene

Inner Core: Polypropylene

Outer Cage: Polypropylene

Sealing Method: Thermal Bonding

O-rings: Buna, Viton(or FKM), EP, Silicone,
Encapsulated Silicone,
Encapsulated Viton (or FKM)

Food Safety Compliance

Material of construction comply with FDA regulations for food and beverage contact use as detailed in the US Code of Federal Regulation, 21CFR.

Material used to produce filter media and hardware are safe for use in contact with foodstuffs in accordance with EU Directives 10/2011.

Eg.=> **CMKPS0201SHHOSW6PBY0**

ORDERING INFORMATION										
Product Type	Membrane Type	Removal Rating	Application	Size	Diameter	Endcap	Inner Core	Sealing Material	Connection Support	Revision
CMK= Mini Pleated Cartridge	PS = PES	0201=0.1/0.2um	S = Ster Grade	HH=500cm ²	OS=STD	W6=116 w/o locking groove Y6=116 w/ locking groove 08=118 double o-ring W8=118 w/ short lugs Y8=118 w/longlugs T2=1/2"NPT threaded 03=123 double o-ring 06=126 double o-ring 1TC=1"TC M8= 80mm skirting T1=1/2"TC B1=1/2"BSP-T	P=PP	B=Buna-N	Y= With SS reinforced	0 = Rev.0
		0601=0.1/0.65um		SS=1000cm ²				E=EPDM		
		0402=0.2/0.45um		DD=1500cm ²				S=Silicone		
		0602=0.2/0.65um		LL=2100cm ²				V=Viton		
		0604=0.45/0.65um								
		0804=0.45/0.8um								

Bio-burden Reduction Junior Cartridge filters

Hydrophobic PTFE membrane

CMK Grade Filter hydrophobic PTFE cartridge are 100% integrity tested; it is constructed of single layer sterilizing expanded polytetrafluoroethylene membrane. It offers broad chemical compatibility, higher filter area with high flow rates at low pressure drops, and low extractables.



Benefits

- 100% integrity tested
- FDA food contact compliant
- Thermal bonding
- Non-fiber releasing

Filter Area

HH: 500cm² SS: 1000cm²
 DD: 1500cm² LL: 2100cm²

Integrity Test Specifications

Pore Size	0.2um
Bubble Point	≥1.4 barg (IPA/Water)
Water Intrusion (area:2100cm ²)	≤0.17ml/min @ 2500mbar, 25°C

Sanitization/Sterilization

- In-line Steam 135°C(275°F), 30 min, 150 cycles

Typical Applications

- Fementer
- Formulation Tank air feed
- RO and WFI storage tank
- Fermentation purge air
- Container purge air
- Aspetic packaging air

Maximum Operating Parameters

- Differential Pressure
 - Forward 50 psid (3.4 bard) at 20°C (68°F)
 - Reverse 40 psid (2.7 bard) at 20°C (68°F)
- Operating Temperature 82°C(180°F) at 10 psid (0.69 barg) in water
- Recommended Changeout Pressure: 35psid(2.4barg)

Toxicity

All components meet the specifications for biological safety per USP Class VI-121°C for plastics.

Construction Materials

- Filter Membrane: PTFE Membrane
- Media Support: Polypropylene
- End Caps: Polypropylene with embedded 316 ring
- Inner Core: SUS 316L perforated core
- Outer Cage: Polypropylene
- Sealing Method: Thermal Bonding
- O-rings: Buna, Viton(or FKM), EP, Silicone,
Encapsulated Silicone,
Encapsulated Viton (or FKM)

Food Safety Compliance

Material of construction comply with FDA regulations for food and beverage contact use as detailed in the US Code of Federal Regulation, 21CFR.
Material used to produce filter media and hardware are safe for use in contact with foodstuffs in accordance with EU Directives 10/2011.

Eg.=> **CMKPT0020BHHOSW6PBY0**

ORDERING INFORMATION										
Product Type	Membrane Type	Removal Rating	Application	Size	Diameter	Endcap	Inner Core	Sealing Material	Connection Support	Revision
CMK=	PT =	0020=0.2um	B =	HH=500cm ²	OS=STD	W6=116 w/o locking groove	P=PP	B=Buna-N	Y= With SS	0 = Rev.0
Mini	PTFE		Low	SS=1000cm ²		Y6=116 w/ locking groove		E=EPDM	reinforced	
Pleated	phobic		Bio	DD=1500cm ²		08=118 double o-ring		S=Silicone		
Cartridge				LL=2100cm ²		W8=118 w/ short lugs		V=Viton		
						Y8=118 w/longlugs				
						T2=1/2"NPT threaded				
						03=123 double o-ring				
						06=126 double o-ring				
						1TC=1"TC				
						M8= 80mm skirting				
						T1=1/2"TC				
						B1=1/2"BSP-T				

General Application Junior Cartridge filters

Hydrophobic PTFE membrane

CMK Grade Filter hydrophobic PTFE cartridge are 100% integrity tested; it is constructed of single layer sterilizing expanded polytetrafluoroethylene membrane. It offers broad chemical compatibility, higher filter area with high flow rates at low pressure drops, and low extractables.



Benefits

- 100% integrity tested
- FDA food contact compliant
- Thermal bonding
- Non-fiber releasing

Typical Applications

- Compressed air
- Pressurized gases
- Tank ventilation
- Chemicals
- Solvents

Filter Area

HH: 500cm² SS: 1000cm²
 DD: 1500cm² LL: 2100cm²

Integrity Test Specifications

Pore size	Bubble point/IPA
0.1um	1.7barg
0.2um	1.1barg
0.45um	0.6barg
1.0um	0.4barg
3.0um	0.1barg
5.0um	0.07barg

Maximum Operating Parameters

- Differential Pressure
 - Forward 80 psid (5.5 barg) at 20°C (68°F)
 - Reverse 40 psid (3.1 barg) at 44°C (75°F)
- Operating Temperature 82°C(180°F) at 10 psid (6.9 barg)
- Recommended Changeout Pressure: 35psid(2.4barg)

Sanitization/Sterilization

- In-line Steam 135°C(275°F), 30 min, 150 cycles

Toxicity

All components meet the specifications for biological safety per USP Class VI-121°C for plastics.

Construction Materials

- Filter Membrane: PTFE Membrane
- Media Support: Polypropylene
- End Caps: Polypropylene with embedded 316 ring
- Inner Core: SUS 316L perforated core
- Outer Cage: Polypropylene
- Sealing Method: Thermal Bonding
- O-rings: Buna, Viton(or FKM), EP, Silicone,
Encapsulated Silicone,
Encapsulated Viton (or FKM)

Food Safety Compliance

Material of construction comply with FDA regulations for food and beverage contact use as detailed in the US Code of Federal Regulation, 21CFR.
Material used to produce filter media and hardware are safe for use in contact with foodstuffs in accordance with EU Directives 10/2011.

Eg.=> **CMKPT0010GHHOSW6PBY0**

ORDERING INFORMATION										
Product Type	Membrane Type	Removal Rating	Application	Size	Diameter	Endcap	Inner Core	Sealing Material	Connection Support	Revision
CMK= Mini Pleated Cartridge	PT = PTFE phobic	0010=0.1um	G = Gen Purpose	HH=500cm ²	OS=STD	W6=116 w/o locking groove Y6=116 w/ locking groove 08=118 double o-ring W8=118 w/ short lugs Y8=118 w/longlugs T2=1/2"NPT threaded 03=123 double o-ring 06=126 double o-ring 1TC=1"TC M8= 80mm skirting T1=1/2"TC B1=1/2"BSP-T	P=PP	B=Buna-N	Y= With SS reinforced	0 = Rev.0
		0020=0.2um		SS=1000cm ²				E=EPDM		
		0045=0.45um		DD=1500cm ²				S=Silicone		
		0065=0.65um		LL=2100cm ²				V=Viton		
		0100=1.0um								
		0300=3.0um								
		0500=5.0um								

General Application Junior Cartridge filters

Nylon 6,6 membrane

CMK Pleated Cartridge Filters with Nylon 66 membrane are 100% integrity tested; it is constructed of single layer hydrophilic nylon66 membrane. It offers broad chemical compatibility, higher filter area with high flow rates at low pressure drops, and low extractables.



Benefits

- 100% integrity tested
- FDA food contact compliant
- Thermal bonding
- Non-fiber releasing

Typical Applications

- Bottled water
- Vinegar
- Aqueous Solutions
- Beer, Wine and Spirits
- Juice, Soft Drinks, Edible Oils
- Pharmaceutical intermediates

Integrity Test Specifications

Wetted membrane

Pore size	Minimum Bubble point
0.1um	50.8 psig (3.5 bar)
0.2um	40.6 psig (2.8 bar)
0.45um	26.1 psig (1.8 bar)
1.0um	8.7 psig (0.6 bar)

Sanitization/Sterilization

- Filtered Hot water 90°C(194°F), 30 min, multiple cycles, max 3 psid forward flow
- Autoclave 121°C(250°F), 30 min, multiple cycles

Maximum Operating Parameters

- Differential Pressure
 - Forward 80 psid (5.5 barg) at 20°C (68°F)
 - Reverse 40 psid (3.1 barg) at 44°C (75°F)
- Operating Temperature 82°C(180°F) at 10 psid (6.9 barg)
- Recommended Changeout Pressure: 35psid(2.4barg)

Toxicity

All components meet the specifications for biological safety per USP Class VI-121°C for plastics.

Construction Materials

Filter Membrane: Nylon Membrane
 Media Support: Polypropylene
 End Caps: Polypropylene with embedded 316 ring
 Inner Core: SUS 316L perforated core
 Outer Cage: Polypropylene
 Sealing Method: Thermal Bonding
 O-rings: Buna, Viton(or FKM), EP, Silicone,
 Encapsulated Silicone,
 Encapsulated Viton (or FKM)

Food Safety Compliance

Material of construction comply with FDA regulations for food and beverage contact use as detailed in the US Code of Federal Regulation, 21CFR.
 Material used to produce filter media and hardware are safe for use in contact with foodstuffs in accordance with EU Directives 10/2011.

Eg.=> **CMKNY0010GHHOSW6PBY0**

ORDERING INFORMATION										
Product Type	Membrane Type	Removal Rating	Application	Size	Diameter	Endcap	Inner Core	Sealing Material	Connection Support	Revision
CMK= Mini Pleated Cartridge	NY = Nylon	0010=0.1um	G = Gen Purpose	HH=500cm ²	OS=STD	W6=116 w/o locking groove Y6=116 w/ locking groove 08=118 double o-ring W8=118 w/ short lugs Y8=118 w/longlugs T2=1/2"NPT threaded 03=123 double o-ring 06=126 double o-ring 1TC=1"TC M8= 80mm skirting T1=1/2"TC B1=1/2"BSP-T	P=PP	B=Buna-N	Y= With SS reinforced	0 = Rev.0
		0020=0.2um		SS=1000cm ²				E=EPDM		
		0045=0.45um		DD=1500cm ²				S=Silicone		
		0065=0.65um		LL=2100cm ²				V=Viton		
		0100=1.0um								

General Application Junior Cartridge filters

Hydrophilic PTFE membrane

CMK Pleated Cartridge Filters with hydrophilic PTFE membrane are 100% integrity tested; it is constructed of single layer hydrophilic PTFE membrane. It offers broad chemical compatibility, higher filter area with high flow rates at low pressure drops, and low extractables.



Benefits

- 100% integrity tested
- FDA food contact compliant
- Thermal bonding
- Non-fiber releasing

Typical Applications

- Bottled water
- Vinegar
- Aqueous Solutions
- Beer, Wine and Spirits
- Juice, Soft Drinks, Edible Oils
- Pharmaceutical intermediates

Integrity Test Specifications

Wetted membrane

Pore size	Minimum Bubble point
0.1um	65.3 psig (4.5 bar)
0.2um	46.4 psig (3.2 bar)
0.45um	31.9 psig (2.2 bar)
1.0um	8.7 psig (0.6 bar)

Sanitization/Sterilization

- Filtered Hot water 90°C(194°F), 30 min, multiple cycles, max 3 psid forward flow
- Autoclave 121°C(250°F), 30 min, multiple cycles

Maximum Operating Parameters

- Differential Pressure
 - Forward 80 psid (5.5 barg) at 20°C (68°F)
 - Reverse 40 psid (3.1 barg) at 44°C (75°F)
- Operating Temperature 82°C(180°F) at 10 psid (6.9 barg)
- Recommended Changeout Pressure: 35psid(2.4barg)

Toxicity

All components meet the specifications for biological safety per USP Class VI-121°C for plastics.

Construction Materials

Filter Membrane: PTFE Membrane
 Media Support: Polypropylene
 End Caps: Polypropylene
 Inner Core: PP/SUS 316L perforated core
 Outer Cage: Polypropylene
 Sealing Method: Thermal Bonding
 O-rings: Buna, Viton(or FKM), EP, Silicone,
 Encapsulated Silicone,
 Encapsulated Viton (or FKM)

Food Safety Compliance

Material of construction comply with FDA regulations for food and beverage contact use as detailed in the US Code of Federal Regulation, 21CFR.
 Material used to produce filter media and hardware are safe for use in contact with foodstuffs in accordance with EU Directives 10/2011.

Eg.=> CMKPH0010GHHOSW6PBYO

ORDERING INFORMATION										
Product Type	Membrane Type	Removal Rating	Application	Size	Diameter	Endcap	Inner Core	Sealing Material	Connection Support	Revision
CMK= Mini Pleated Cartridge	PH = PTFE phillic	0010=0.1um	G = Gen Purpose	HH=500cm ²	OS=STD	W6=116 w/o locking groove Y6=116 w/ locking groove 08=118 double o-ring W8=118 w/ short lugs Y8=118 w/longlugs T2=1/2"NPT threaded 03=123 double o-ring 06=126 double o-ring 1TC=1"TC M8= 80mm skirting T1=1/2"TC B1=1/2"BSP-T	P=PP	B=Buna-N	Y= With SS reinforced	0 = Rev.0
		0020=0.2um		SS=1000cm ²				E=EPDM		
		0045=0.45um		DD=1500cm ²				S=Silicone		
		0065=0.65um		LL=2100cm ²				V=Viton		
		0100=1.0um								

Junior Cartridge filters

Polypropylene media

CMK Cartridge Filters with depth structure of polypropylene media. It offers broad chemical compatibility, higher dirt holding capacity with high flow rates at low pressure drop and low extractables. They are available in nominal and absolute ratings.



Benefits

- Wide chemical compatibility
- High dirt hold capacity
- High retention
- Thermal bonding
- Non-fiber releasing

Typical Applications

- Process water
- Vinegar
- Aqueous Solutions
- Beer, Wine and Spirits
- Juice, Soft Drinks, Edible Oils
- Bulk Chemicals
- Pharmaceutical intermediates

Maximum Operating Parameters

- Differential Pressure
 - Forward 80 psid (5.5 barg) at 20°C (68°F)
 - Reverse 40 psid (3.1 barg) at 44°C (75°F)
- Operating Temperature 82°C(180°F) at 10 psid (6.9 barg)
- Recommended Changeout Pressure: 35psid(2.4barg)

Sanitization/Sterilization

- Filtered Hot water 90°C(194°F), 30 min, multiple cycles, max 3 psid forward flow
- Autoclave 121°C(250°F), 30 min, multiple cycles

Toxicity

All components meet the specifications for biological safety per USP Class VI-121°C for plastics.

Construction Materials

Filter Membrane: PP Membrane
 Media Support: Polypropylene
 End Caps: Polypropylene
 Inner Core: SUS 316L perforated core
 Outer Cage: Polypropylene
 Sealing Method: Thermal Bonding
 O-rings: Buna, Viton(or FKM), EP, Silicone,
 Encapsulated Silicone,
 Encapsulated Viton (or FKM)


Food Safety Compliance

Material of construction comply with FDA regulations for food and beverage contact use as detailed in the US Code of Federal Regulation, 21CFR.
 Material used to produce filter media and hardware are safe for use in contact with foodstuffs in accordance with EU Directives 10/2011.

Eg.=> CMKPP0060GHHOSW6PBY0

ORDERING INFORMATION										
Product Type	Membrane Type	Removal Rating	Application	Size	Diameter	Endcap	Inner Core	Sealing Material	Connection Support	Revision
CMK= Mini Pleated Cartridge	PP = Polypro	0060=0.6um	G = Gen Purpose	HH=500cm ²	OS=STD	W6=116 w/o locking groove Y6=116 w/ locking groove O8=118 double o-ring W8=118 w/ short lugs Y8=118 w/longlugs T2=1/2"NPT threaded O3=123 double o-ring O6=126 double o-ring 1TC=1"TC M8= 80mm skirting T1=1/2"TC B1=1/2"BSP-T	P=PP	B=Buna-N	Y= With SS reinforced	0 = Rev.0
		0100=1.0um		SS=1000cm ²						
		0300=3.0um	P = Premier	DD=1500cm ²				S=Silicone		
		0500=5.0um		LL=2100cm ²				V=Viton		
		0700=7.0um								
		1000=10.0um								
		2000=20.0um								
		3000=30.0um								
5000=50.0um										



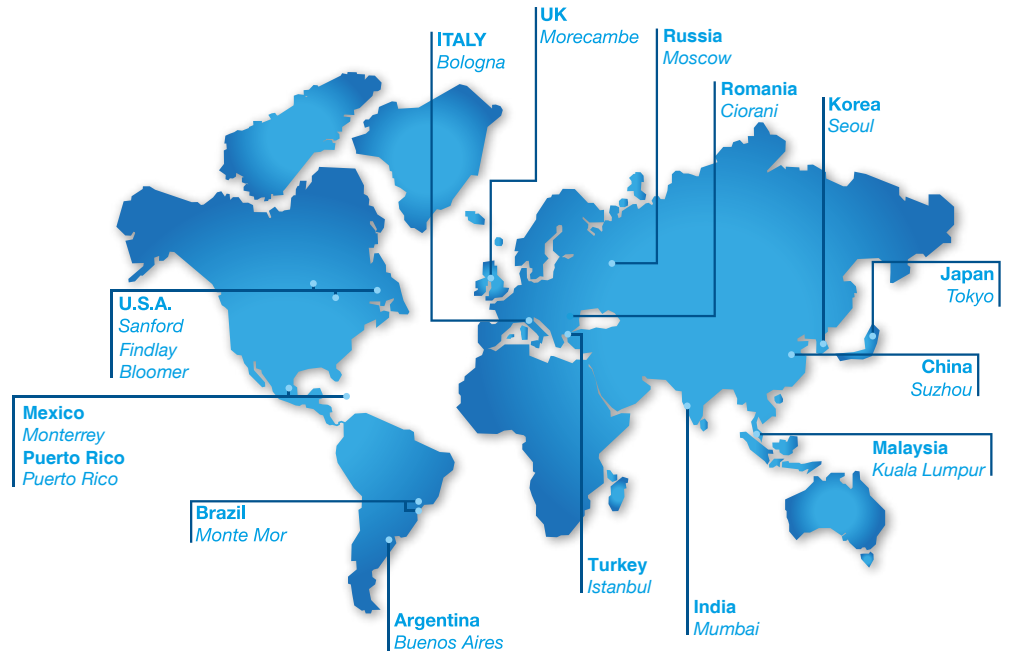

WORLDWIDE

EUROPE

Italy Office
Headquarters
GVS S.p.A.
Via Roma 50
40069 Zola Predosa (BO) - Italy
Tel. +39 051 6176311
gvs@gvs.com

Russia
GVS Russia LLC.
Profsoyuznaya Street, 25-A, office 102
117418, Moscow
Russian Federation (Russia)
Tel. +7 495 0045077
gvsrussia@gvs.com

United Kingdom
GVS Filter Technology UK
Vickers Industrial Estate
Mellishaw Lane, Morecambe
Lancashire LA3 3EN
Tel. +44 (0) 1524 847600
gvsuk@gvs.com

Romania
GVS Microfiltrazione srl
Sat Ciorani de Sus 1E - Comuna Ciorani
Prahova România
Tel. (+40) 244 463044
gvsro@gvs.com

Turkey
GVS Türkiye
Nidakule Merdivenköy Mahallesi
Bora Sokak No:1 Kat:7 - 34732 Istanbul
Tel. +90 216 504 47 67
gvsurkey@gvs.com

ASIA

China
GVS Technology (Suzhou) Co., Ltd.
Fengqiao Civil-Run Sci-Tech Park,
602 Changjiang Road,S.N.D.
Suzhou, China 215129
Tel. +86 512 6661 9880
gvschina@gvs.com

GVS YIBO Medical Devices Co. Ltd.
17, Zhongshan East - Yuyao city,
315403 Zhejiang Province, China
Tel. +86 574 6257 5697

Japan
GVS Japan K.K.
KKD Building 4F, 7-10-12 Nishishinjuku
Shinjuku-ku, Tokyo 160-0023 Japan
Tel. +81 3 5937 1447
gvsjapan@gvs.com

Korea
GVS Korea Ltd
#315 Bricks Tower
368 Gyungchun-ro(Gaun-dong),
Namyangju-si, Gyunggi-do,
Tel: +82 31 563 9873
gvsukorea@gvs.com

India
GVS Filter India Pvt Ltd
Unit No 35 & 36 on First Floor
Ratna Jyot Industrial Premises Irla Lane,
Irla Vile Parle, Mumbai 400056, India

Malaysia
GVS Filtration Sdn.Bhd
Lot No 10F-2B, 10th Floor, Tower 5 @ PFCC
Jalan Puteri 1/2, Bandar Puteri
47100 Puchong, Selangor, Malaysia

AMERICA

U.S.A.
GVS North America
63 Community Drive
Sanford, ME 04073 - USA
Tel. +1 866 7361250
gvsusa@gvs.com

GVS Filtration Inc.
2150 Industrial Drive
Findlay, OH. 45840 - USA
Tel. +1.419.423.9040
gvsfiltration@gvs.com

2200 W 20th Avenue
Bloomer, WI 54724 - USA
Tel. +1.715.568.5944
gvsfiltration@gvs.com

Puerto Rico
GVS Puerto Rico, LLC
98 Carr 194 - Fajardo,
Puerto Rico, 00738-2988, USA
Tel. +1.787.355.4100
gvspuertorico@gvs.com

México
GVS Filter Technology de Mexico
Universal No. 550, Ynmsa Aeropuerto Apodaca
Industrial Park, Ciudad Apodaca, Nuevo León,
C.P. 66626 - México
Tel. +52 81 2282 9003
gvmex@gvs.com

Argentina
GVS Argentina S.A.
Francisco Acuña de Figueroa 719 Piso:11 Of: 57
1416 Buenos Aires - Argentina
Tel. +54 11 48614750
gvsarg@gvs.com

Brazil
GVS do Brasil Ltda.
Rodovia Conego Cyriaco Scaranello Pires 251
Jardim Chapadão, CEP 13193-580
Monte Mor (SP) - Brasil
Tel. +55 19 38797200
gvs@gvs.com.br

PRODUCT COLLECTION - Junior cartridge catalog

Copyright © 2021 GVS ® S.p.A.
All Right Reserved
Printing History: Version 29062022